

Analyte Specifications:

**CD154  
(CD40L)**

Catalog Number:


**92008**

**MultiCyt™ QBeads™ PlexScreen**

Reconstituting the CD154 (CD40L) standard in 2 mL of assay diluent results in a concentration of 5,000 pg/mL. Standard curves for no-wash and 1-wash protocols were performed in triplicate (+/- S.D.). The limit of detection (dashed line) is 3 standard deviations above the mean of samples without Standard (N=12). Results are typical and may vary based on specific assay conditions.


### STANDARD CURVES

**No-Wash Assay Performance**

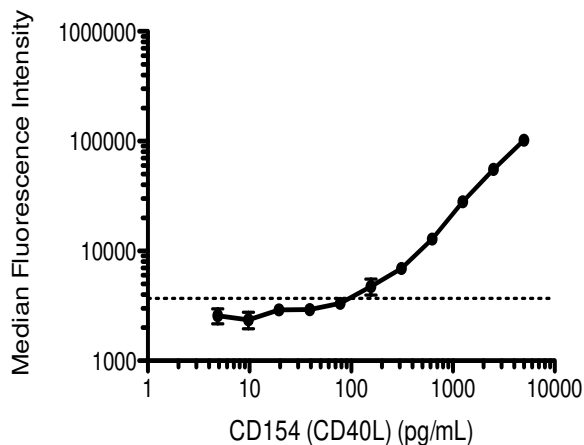

**Lowest Detectable Test Conc. (pg/mL): 156.3**

**1-Wash Assay Performance**

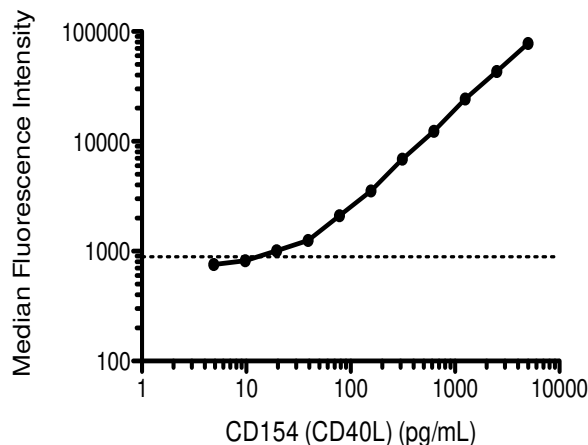

**Lowest Detectable Test Conc. (pg/mL): 19.5**

For research purposes only. Not for resale.

Rev. 1 Jan 2014