

Analyte Specifications:  
**CD62P**  
**(P-Selectin)**


Catalog Number:  
**92006**

**MultiCyt™ QBeads™ PlexScreen**


Reconstituting the CD62P standard in 2 mL of assay diluent results in a concentration of 20,000 pg/mL. Standard curves for no-wash and 1-wash protocols were performed in triplicate (+/- S.D.). The limit of detection (dashed line) is 3 standard deviations above the mean of samples without Standard (N=12). Results are typical and may vary based on specific assay conditions.

## STANDARD CURVES

**No Wash Protocol**

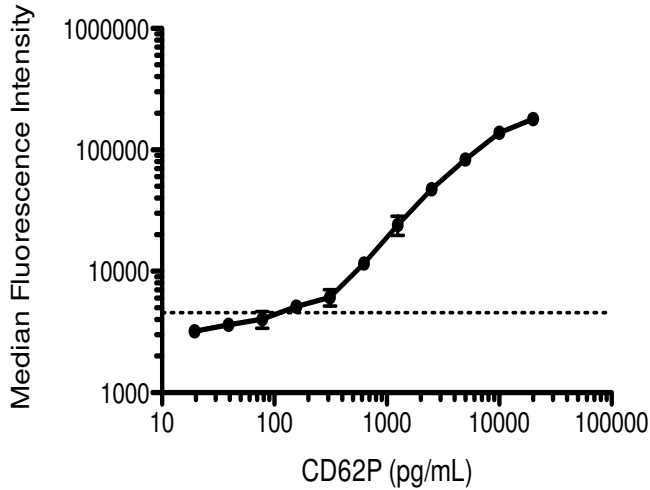

**1-Wash Protocol**

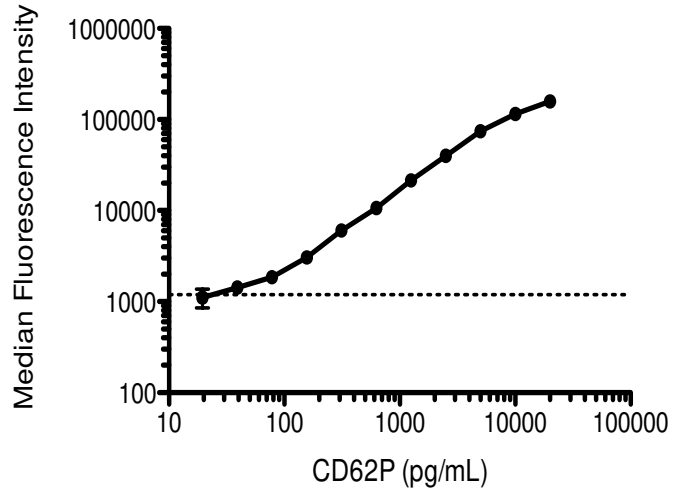

**Lowest Detectable Test Conc. (pg/mL): 156.3**

**Lowest Detectable Test Conc. (pg/mL): 39.1**

For research purposes only. Not for resale.

Rev. 1 Jan 2014