

Analyte Specifications:  
**CX3CL1**  
**(Fractalkine)**


Catalog Number:  
**92011**

**MultiCyt™ QBeads™ PlexScreen**

Reconstituting the Fractalkine standard in 2 mL of assay diluent results in a concentration of 20,000 pg/mL. Standard curves for no-wash and wash protocols were performed in triplicate (+/- S.D.). The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without Fractalkine (N=12). Results are typical and may vary based on specific assay conditions.


## STANDARD CURVES

**No Wash Protocol**

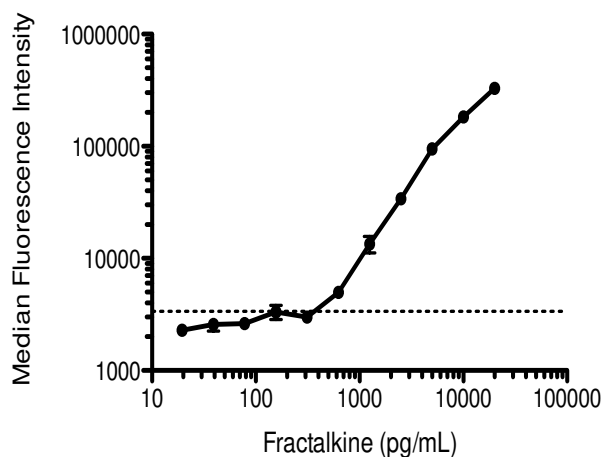

**Lowest Detectable Test Conc. (pg/mL): 625**

**1-Wash Protocol**

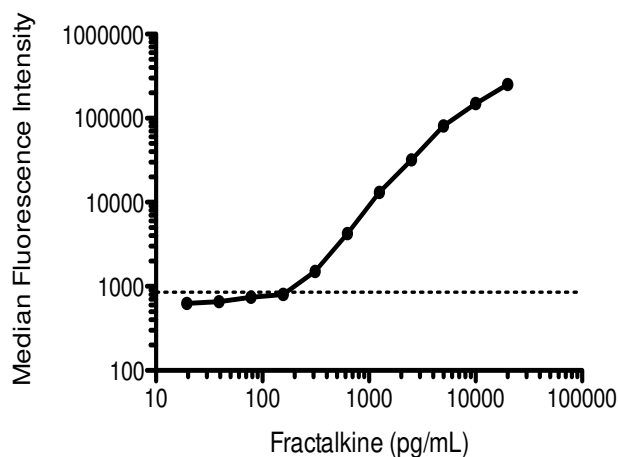

**Lowest Detectable Test Conc. (pg/mL): 312.5**

For research purposes only. Not for resale.

Rev. 1 Jan 2014