

Analyte Specifications:
CSF2
(GM-CSF)


Catalog Number:
92014

MultiCyt™ QBeads™ PlexScreen

Reconstituting the GM-CSF standard in 2 mL of assay diluent results in a concentration of 5000 pg/mL. Standard curves for no-wash and wash protocols were run in triplicate (+/- S.D.). The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without GM-CSF (N=12). Results are typical and may vary based on specific assay conditions.


STANDARD CURVES

No Wash Protocol

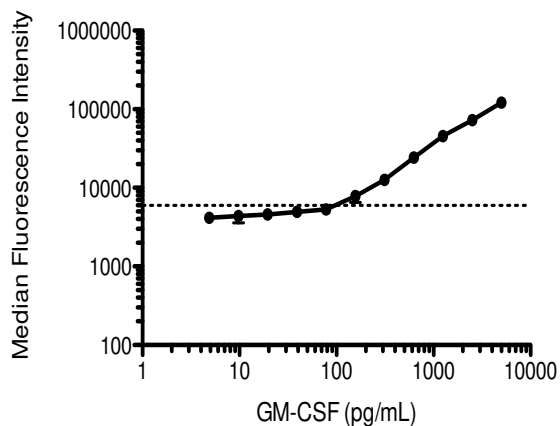

Lowest Detectable Test Conc. (pg/mL): 156.3

1-Wash Protocol

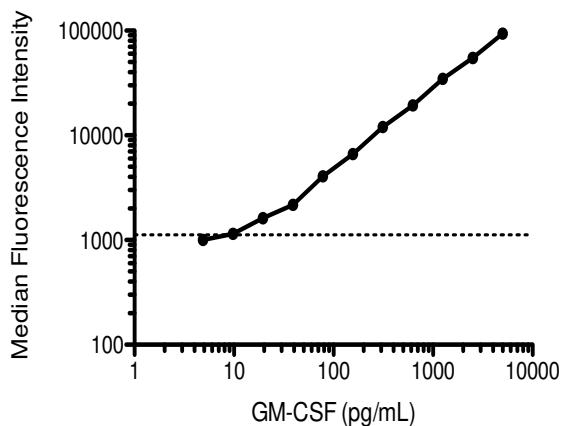

Lowest Detectable Test Conc. (pg/mL): 9.8

For research purposes only. Not for resale.

Rev. 1 Jan 2014