

Analyte Specifications:

Catalog Number:

**IL-12/IL-23p40**


**92029**

**MultiCyt™ QBeads™ PlexScreen**


Reconstituting the IL-12/IL-23p40 standard in 2 mL of assay diluent results in a concentration of 20,000 pg/mL. Standard curves for no-wash and wash protocols were performed in triplicate (+/- S.D.). The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without IL-12/IL-23p40 (N=12). Results are typical and may vary based on specific assay conditions.

## STANDARD CURVES

**No Wash Protocol**

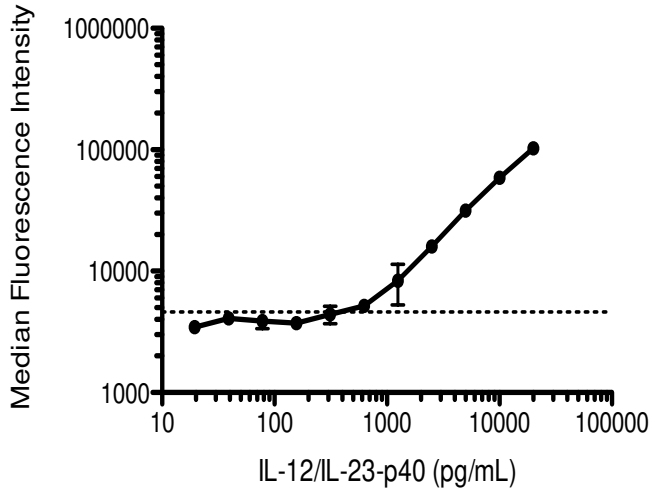

**1-Wash Protocol**

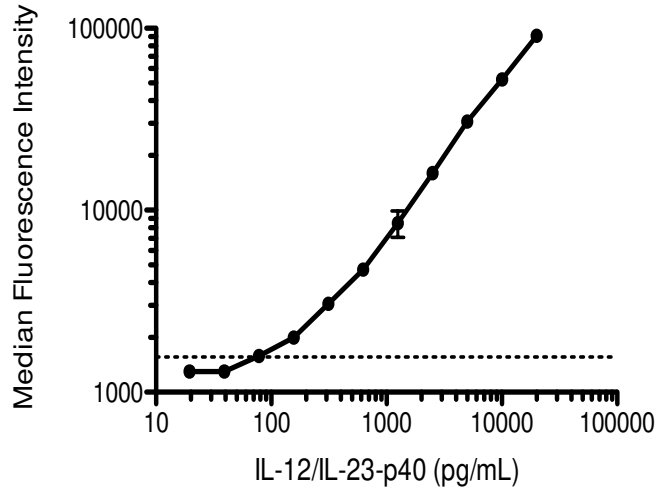

**Lowest Detectable Test Conc. (pg/mL): 625**

**Lowest Detectable Test Conc. (pg/mL): 78.1**

For research purposes only. Not for resale.