

Analyte Specifications:

Catalog Number:

IL-4


92021

MultiCyt™ QBeads™ PlexScreen


Reconstituting the IL-4 standard in 2 mL of assay diluent results in a concentration of 5000 pg/mL. Standard curves for no-wash and wash protocols were performed in triplicate (+/- S.D.) The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without IL-4 (N=12). Results are typical and may vary based on specific assay conditions.

STANDARD CURVES

No Wash Protocol

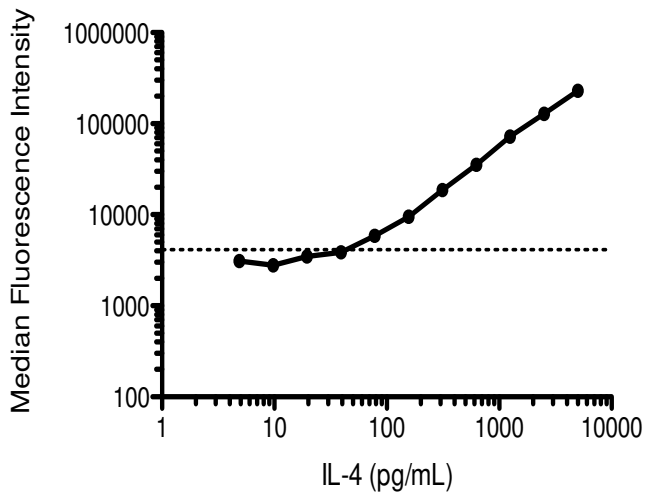

1-Wash Protocol

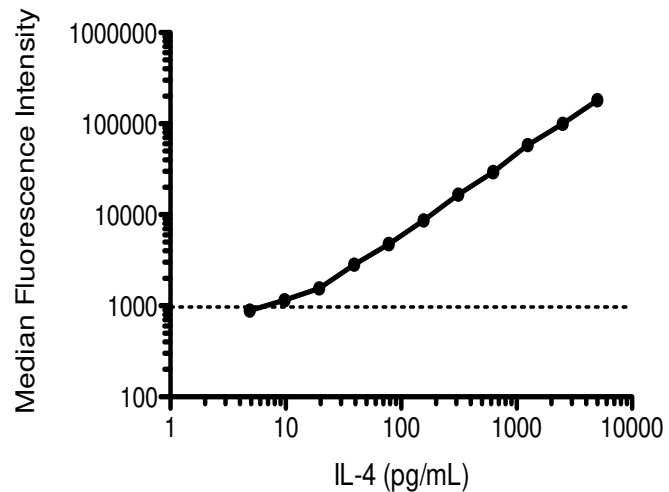

Lowest Detectable Test Conc. (pg/mL): 78.1

Lowest Detectable Test Conc. (pg/mL): 9.8

For research purposes only. Not for resale.