

Analyte Specifications:

Catalog Number:

**CSF3
(G-CSF)**



92054

**MultiCyt® QBeads™
PlexScreen Mouse Secreted
Protein**

Reconstituting the G-CSF standard in 2 mL of assay diluent results in a concentration of 5,000 pg/mL. Standard curves for no-wash and 1-wash protocols are shown starting from 5,000 pg/mL and proceeding in a 1:2 serial titration. The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without G-CSF (N=12). Results are typical and may vary based on specific assay conditions.

No-Wash Assay Performance

1-Wash Assay Performance

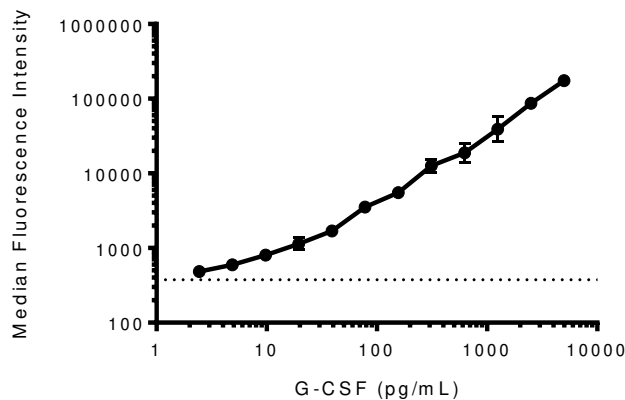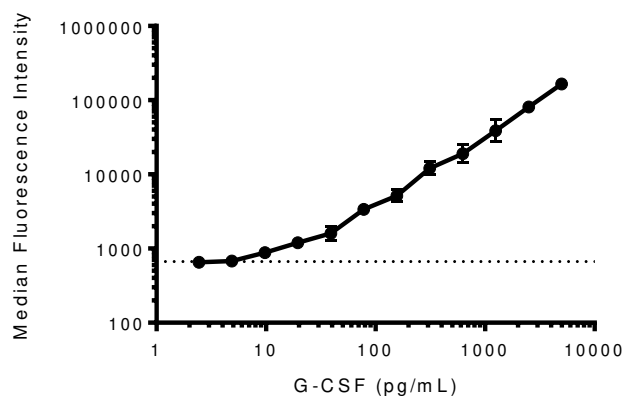

**Lowest Detectable
Concentration Tested (pg/mL): 5**

**Lowest Detectable
Concentration Tested (pg/mL): 2.4**

For research purposes only. Not for resale.

Rev. 1 Sep 2014

