

Analyte Specifications:

Catalog Number:

**IFN $\gamma$**



**92056**

**MultiCyt<sup>®</sup> QBeads<sup>™</sup>  
PlexScreen Mouse Secreted  
Protein**

Reconstituting the IFN $\gamma$  standard in 2 mL of assay diluent results in a concentration of 5,000 pg/mL. Standard curves for no-wash and 1-wash protocols are shown starting from 5,000 pg/mL and proceeding in a 1:2 serial titration. The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without IFN $\gamma$  (N=12). Results are typical and may vary based on specific assay conditions.

**No-Wash Assay Performance**

**1-Wash Assay Performance**

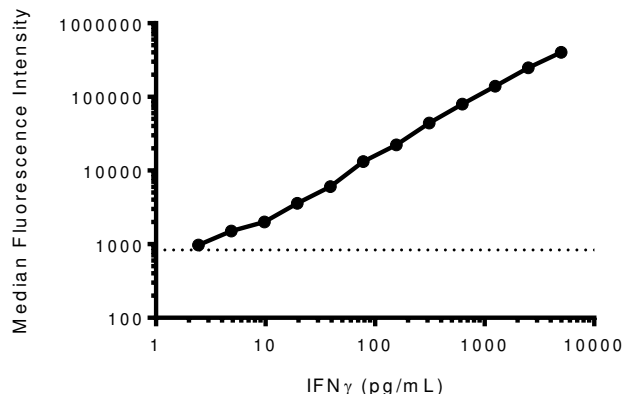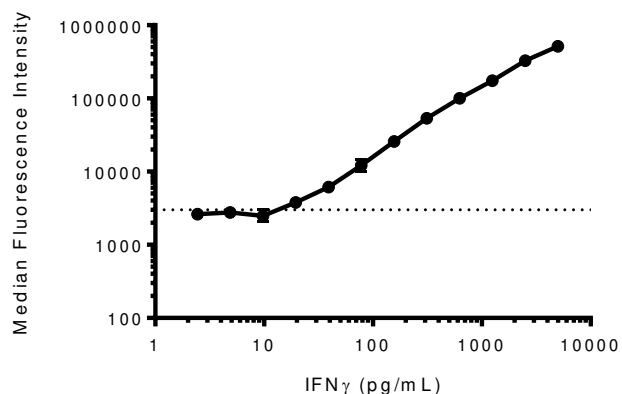

**Lowest Detectable  
Concentration Tested (pg/mL): 20**

**Lowest Detectable  
Concentration Tested (pg/mL): 2.4**

For research purposes only. Not for resale.

Rev. 1 Sep 2014