

Analyte Specifications:

**IL-12/IL-23
(p40)**

Catalog Number:



92081

**MultiCyt® QBeads™
PlexScreen Mouse Secreted
Protein**

Reconstituting the IL-12/IL-23 (p40) standard in 2 mL of assay diluent results in a concentration of 5,000 pg/mL. Standard curves for no-wash and 1-wash protocols are shown starting from 5,000 pg/mL and proceeding in a 1:2 serial titration. The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without IL-12/IL-23 (p40) (N=12). Results are typical and may vary based on specific assay conditions.

No-Wash Assay Performance

1-Wash Assay Performance

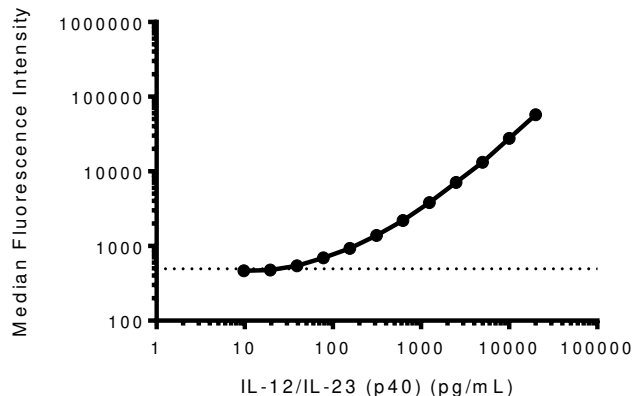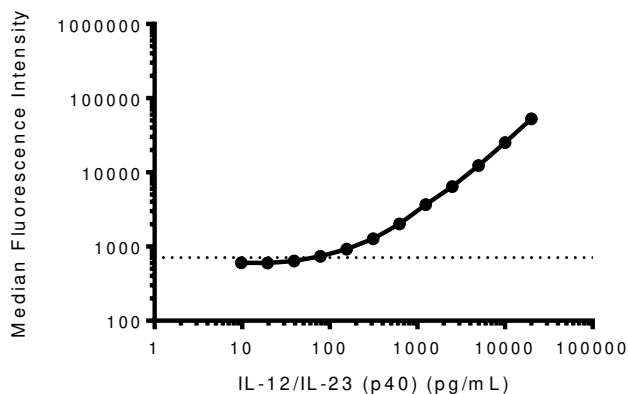

**Lowest Detectable
Concentration Tested (pg/mL): 78**

**Lowest Detectable
Concentration Tested (pg/mL): 39**

For research purposes only. Not for resale.

Rev. 1 Sep 2014

