

Analyte Specifications:

MIP-1 α
(CCL3)

Catalog Number:



92075

MultiCyt[®] QBeads[™]
PlexScreen Mouse Secreted
Protein

Reconstituting the MIP-1 α standard in 2 mL of assay diluent results in a concentration of 5000 pg/mL. Standard curves for no-wash and 1-wash protocols are shown starting from 5000 pg/mL and proceeding in a 1:2 serial titration. The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without MIP-1 α (N=12). Results are typical and may vary based on specific assay conditions.

No-Wash Assay Performance

1-Wash Assay Performance

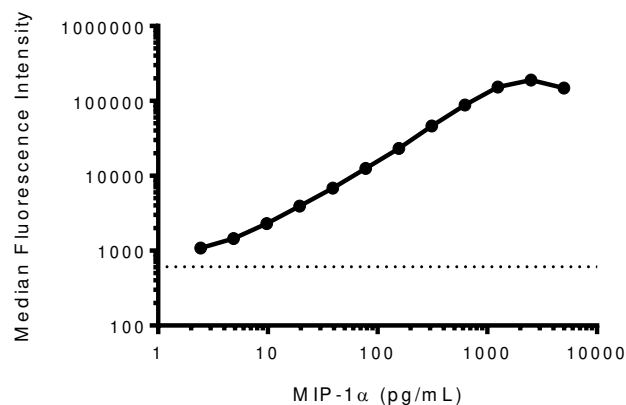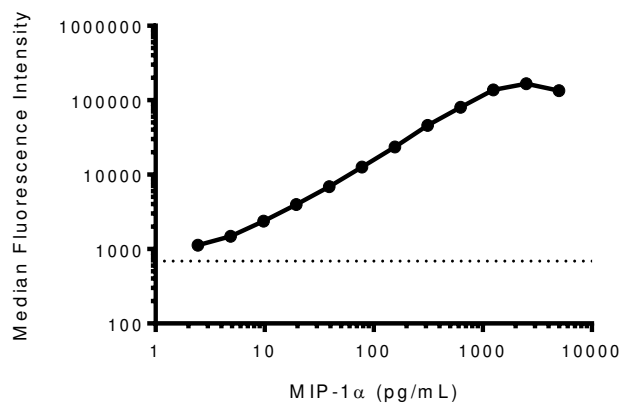

Lowest Detectable

Concentration Tested (pg/mL): 2.4

Lowest Detectable

Concentration Tested (pg/mL): 2.4

For research purposes only. Not for resale.

Rev. 1 Sep 2014