

Analyte Specifications:  
**TNF**  
**(Option 1)**


Catalog Number:  
**92047**

**MultiCyt™ QBeads™ PlexScreen**

Reconstituting the TNF $\alpha$  standard in 2 mL of assay diluent results in a concentration of 5000 pg/mL. Standard curves for no-wash and wash protocols were run in triplicate ( $\pm$  S.D.). The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without TNF $\alpha$  (N=12). Results are typical and may vary based on specific assay conditions.

## STANDARD CURVES


**No Wash Protocol**

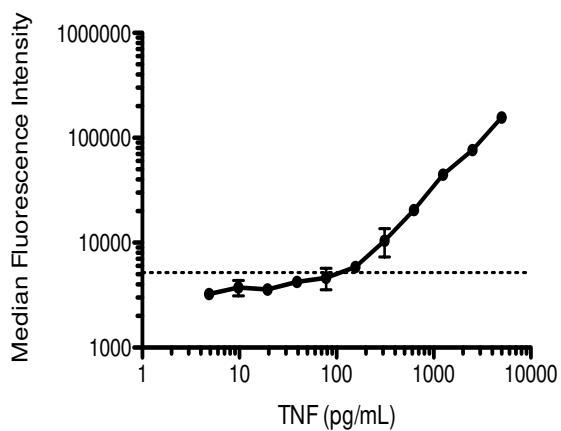

**Lowest Detectable Test Conc. (pg/mL): 156**

For research purposes only. Not for resale.

**1-Wash Protocol**

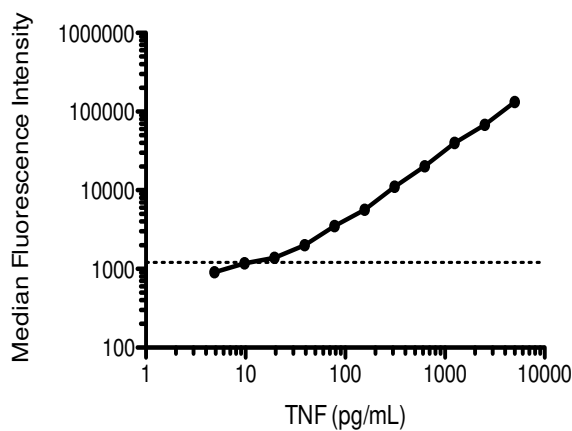

**Lowest Detectable Test Conc. (pg/mL): 20**

Rev 1 Jan. 2014