

Analyte Specifications:

**CD54
(ICAM-1)**

Catalog Number:


92003

MultiCyt™ QBeads™ PlexScreen


Reconstituting the CD54 (ICAM-1) standard in 2 mL of assay diluent results in a concentration of 20,000 pg/mL. Standard curves for no-wash and wash protocols were performed in triplicate (+/- S.D.). The limit of detection, represented by the dashed line, is determined by adding 3 standard deviations to the mean of samples without Standard (N=12). Results are typical and may vary based on specific assay conditions.

STANDARD CURVES

No Wash Protocol

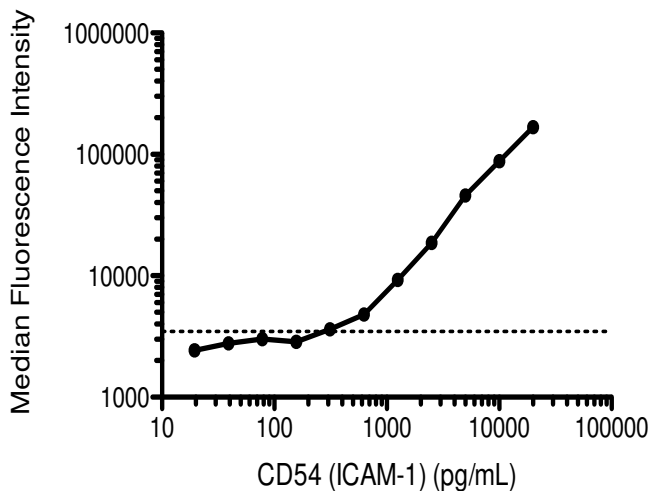

1-Wash Protocol

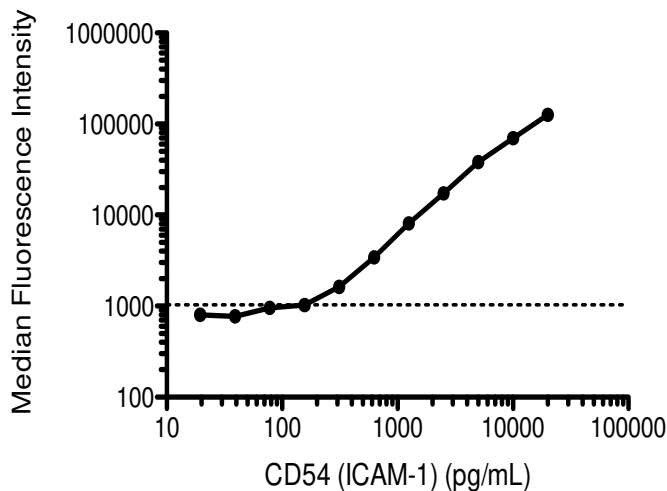

Lowest Detectable Test Conc. (pg/mL): 312

Lowest Detectable Test Conc. (pg/mL): 180

For research purposes only. Not for resale.

Rev. 1 Jan 2014